
GB	Instructions for use	NL	Gebruiksaanwijzing	DE	Gebrauchsanweisung
AR	تعليمات الاستخدام	BA	Uputstvo za upotrebu	CZ	Návod k použití
DK	Brugsanvisning	ES	Instrucciones de uso	FI	Käyttöohje
FR	Mode d'emploi	GR	Οδηγίες χρήσης	HR	Upute za uporabu
ID	Petunjuk Penggunaan	IT	Istruzioni per l'uso	LT	Naudojimo instrukcijos
LV	Lietošanas instrukcija	NO	Bruksanvisning	PL	Instrukcja użycia
PT	Instruções de uso	RO	Instructiuni de folosire	SE	Bruksanvisning
SI	Navodila za uporabo	SK	Inštrukcie na používanie	TR	Kullanım talimatları
UA	Інструкція із застосування				

REF

C1FB108-30	C1FZ108-30	C1NB10	C1NB10-30
C1NB10-HC	C1NB108-30	C1NB25	C1NB29
C1NB32	C1NB38	C1NB10V	C1NB10V-HC-30
C1NB108V-HC-30	C1NB25V	C1NB29V	C1NB32V
C1NB38V	C1NC10	C1NZ10-30	C1NZ108-30
C1NZ25	C1NZ29	C1NZ32	C1NZ38
C1NZ10V	C1NZ10V-HC-30	C1NZ108V-HC-30	C1NZ25V
C1NZ29V	C1NZ32V	C1NZ38V	C1PB10
C1PB10-30	C1PB10V-30	C1PZ10-30	C1PZ10V-30
N1NB10	N1NZ10		

MD 0°C 27°C CE

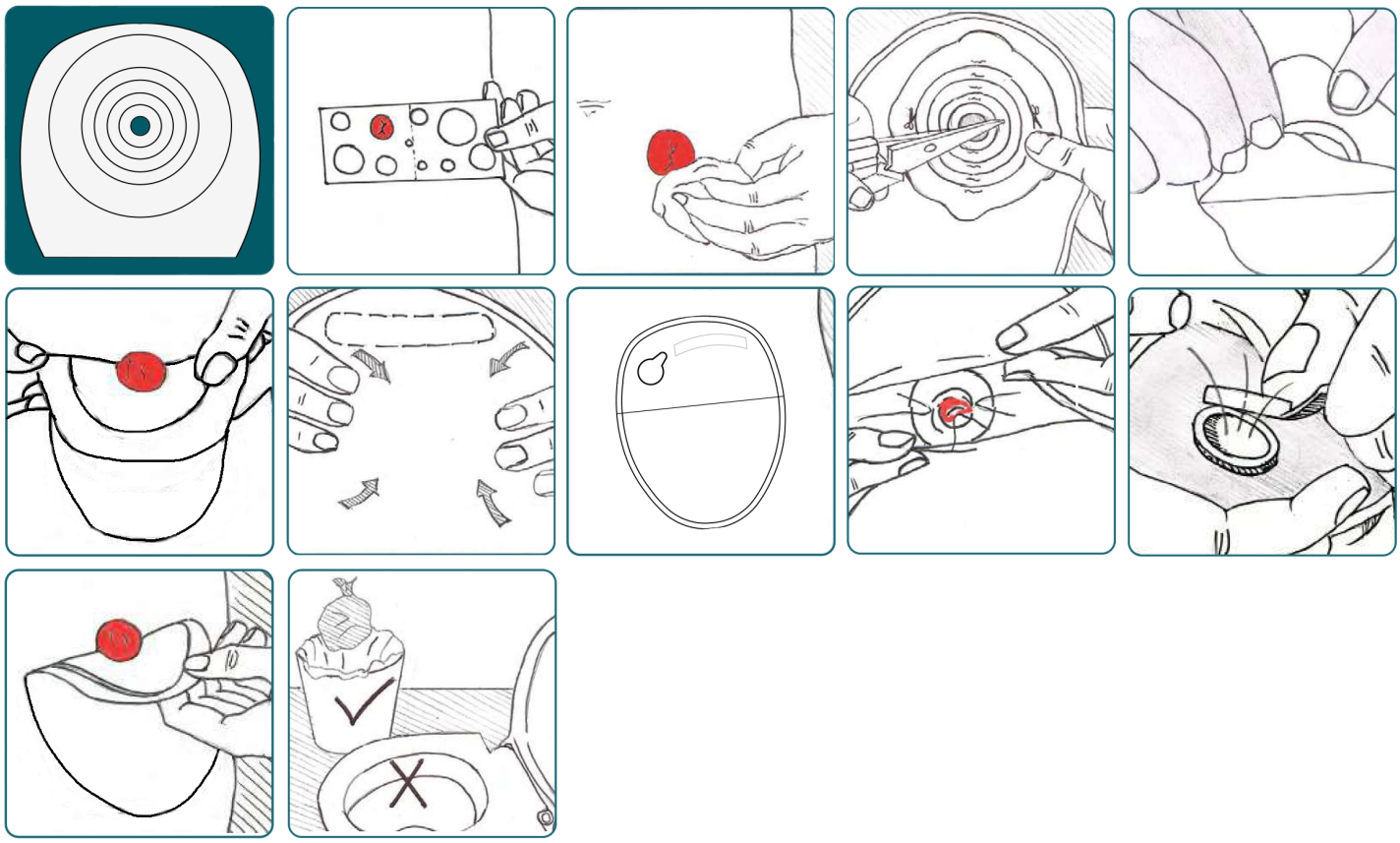GB - Instructions for use

Intended Use, Product Description and Material Information
Ostomy pouches are used to collect waste from the body and are usually attached to the skin over a surgically made stoma. Zensiv Ostomy products are a range of 1-piece and 2-piece ostomy pouches and 2-piece flanges in various sizes and shapes. All 1-piece pouches and 2-piece flanges feature a skin friendly hydrocolloid wafer for adhesion to the body and this is protected before use with a release film that shows a cutting guide. The product pouch is manufactured from odour barrier film with a comfort non-woven layer and a carbon filter for odour control (for ileostomy and colostomy). Ileostomy pouches have either a roll up drain port secured with an integral 'velcro' type closure or there are clamps and soft ties supplied in each box. For urostomy the pouches carry a soft polymer tap for easy drainage. All pouches have a viewing window to inspect the stoma if needed. Some ranges feature a vent with lid for venting the pouch if required.All Zensiv Ostomy products are latex-free and do not contain ingredients of animal origin.The products are non-sterile and can be used in any appropriate circumstances where the attending physician or healthcare professional decides a non-sterile product may be applied.

Warnings-Precautions
- Each pouch is for single use only. Do not re-use. Re-use can lead to cross-contamination. - Use only as directed by your healthcare professional. - Read the instructions prior to use. - Discontinue use and contact your healthcare professional if you experience any allergic reactions or problems when using the pouches.

Contra-indications
These devices should not be used by anyone with a known sensitivity to: polyethylene, polyamide, medical charcoal, synthetic polymer, ethylene vinyl acetate, polypropylene, polyester, polyisobutylene, carboxymethylcellulose, pectin or glycerine, acrylic gel, spandex.

Storage Conditions
It is best to store your products in a cool dry place with a guidance of temperature from 0°C to 27°C as is noted on the product label.

Serious Incidents
Should a serious incident arise as a result of using the product please report this immediately to the manufacturer and the Competent Authority in the country where you reside.

Shelf Life
Use by the expiry date noted on the label.

Collect your materials
Make sure that all necessary materials are ready for use. Collect a measuring guide, a pair of scissors, a pen, wet and dry gauze and a stoma pouch.

Choose the right size
Measure your stoma with a measuring guide and cut (if necessary) the flange to the right size of the stoma. Make sure that you choose the right size so that the stoma and surrounding skin will stay healthy. Allow a maximum of 2 mm of space between your stoma and flange for a safe and correct fit. If the size is too big, your faeces will damage the skin. If the size is too small, the flange will damage the stoma. Over time your stoma may change shape or size, so periodically measurements may be required again to prevent any leakage.

Prepare your skin
Support the skin and remove the flange carefully, do not tear off but push the skin back from the flange, use a wet gauze if necessary. Wash the skin surrounding your stoma with warm and unperfumed water. Do not use anything that contains alcohol, perfume or ether to clean your skin as it can cause damage. Do not rub the skin too hard as this might lead to skin damage. Dry the skin carefully and thoroughly to ensure that your pouch adheres well. Optional: If your skin is irritated and sore you can use skin barrier cream, hydrocolloid powder or skin barrier wipes. Apply the powder, skin barrier cream or wipe gently to the surrounding skin.

Fit your pouch
Remove the release paper of the flange. Fold the flange halfway and place the opening at the underside of the stoma. Stretch the skin of your abdomen with a flat hand and pull the whole flange over the stoma into the right position. Work your way around the stoma so that the hydrocolloid will attach to the skin and make sure that the flange is secure. Place your hand over the flange and make sure that a good seal is created.

Vent
If the pouch is full of gasses (ballooning) or if there is no air in the pouch (a vacuum) the vent can be used to let air in or out of the pouch. Make sure that the pouch is not full of faeces when opening the vent. If the pouch is full than empty the pouch or change it for a fresh one. The vent can be opened by taking off the lid. Carefully close the lid after venting.

Viewing window
The ostomy pouch includes a viewing window which makes it possible to inspect the stoma as well as the output from the stoma. It is also very helpful when applying the pouch to make sure it is located correctly over the stoma. Lift up the non-woven of the viewing window for inspection.

Disposal of your pouch
Remove your pouch and dispose of it in the trash. Do not flush your used pouch down the toilet because this will block your sewage system.

Frequency of pouch change
Change/empty your pouch as advised by your healthcare professional.

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Oxmed International GmbH