
| | | | | | |
|---|---|---|---|---|---|
| GB | Instructions for use | NL | Gebruiksaanwijzing | DE | Gebrauchsanweisung |
| AR | تعليمات الاستخدام | BA | Uputstvo za upotrebu | CZ | Návod k použití |
| DK | Brugsanvisning | ES | Instrucciones de uso | FI | Käyttöohje |
| FR | Mode d'emploi | GR | Οδηγίες χρήσης | HR | Upute za uporabu |
| ID | Petunjuk Penggunaan | IT | Istruzioni per l'uso | LT | Naudojimo instrukcijos |
| LV | Lietošanas instrukcija | NO | Bruksanvisning | PL | Instrukcja użycia |
| PT | Instruções de uso | RO | Instructiuni de folosire | SE | Bruksanvisning |
| SI | Navodila za uporabo | SK | Inštrukcie na používanie | TR | Kullanım talimatları |
| UA | Інструкція із застосування | | | | |

REF

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MD 0°C 27°C CE
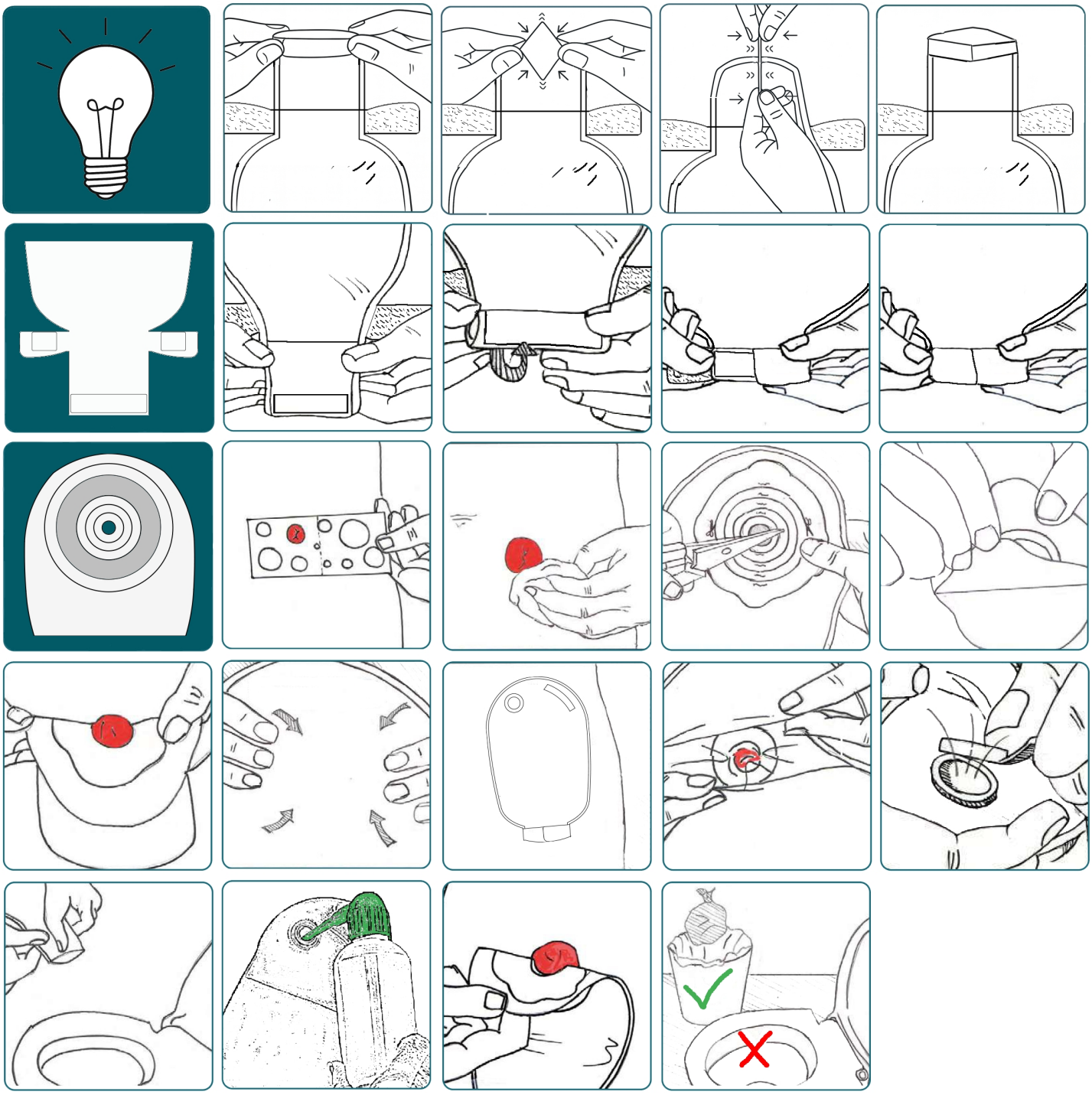
## GB - Instructions for use

**Intended Use, Product Description and Material Information**
Ostomy pouches are used to collect waste from the body and are usually attached to the skin over a surgically made stoma. Zensiv Ostomy products are a range of 1-piece and 2-piece ostomy pouches and 2-piece flanges in various sizes and shapes. All 1-piece pouches and 2-piece flanges feature a skin friendly hydrocolloid wafer for adhesion to the body and this is protected before use with a release film that shows a cutting guide. The product pouch is manufactured from odour barrier film with a comfort non-woven layer and a carbon filter for odour control (for ileostomy and colostomy). Ileostomy pouches have either a roll up drain port secured with an integral 'velcro' type closure or there are clamps and soft ties supplied in each box. For urostomy the pouches carry a soft polymer tap for easy drainage. All pouches have a viewing window to inspect the stoma if needed. Some ranges feature a vent with lid for venting the pouch if required.All Zensiv Ostomy products are latex-free and do not contain ingredients of animal origin.The products are non-sterile and can be used in any appropriate circumstances where the attending physician or healthcare professional decides a non-sterile product may be applied.

**Warnings-Precautions**
- Each pouch is for single use only. Do not re-use. Re-use can lead to cross-contamination. - Use only as directed by your healthcare professional. - Read the instructions prior to use. - Discontinue use and contact your healthcare professional if you experience any allergic reactions or problems when using the pouches.

**Contra-indications**
These devices should not be used by anyone with a known sensitivity to: polyethylene, polyamide, medical charcoal, synthetic polymer, ethylene vinyl acetate, polypropylene, polyester, polyisobutylene, carboxymethylcellulose, pectin or glycerine, acrylic gel, spandex.

**Storage Conditions**
It is best to store your products in a cool dry place with a guidance of temperature from 0°C to 27°C as is noted on the product label.

**Serious Incidents**
Should a serious incident arise as a result of using the product please report this immediately to the manufacturer and the Competent Authority in the country where you reside.

**Shelf Life**
Use by the expiry date noted on the label.

**Collect your materials**
Make sure that all necessary materials are ready for use. Collect a measuring guide, a pair of scissors, a pen, wet and dry gauze and a stoma pouch.

**Choose the right size**
Measure your stoma with a measuring guide and cut (if necessary) the flange to the right size of the stoma. Make sure that you choose the right size so that the stoma and surrounding skin will stay healthy. Allow a maximum of 2 mm of space between your stoma and flange for a safe and correct fit. If the size is too big, your faeces will damage the skin. If the size is too small, the flange will damage the stoma. Over time your stoma may change shape or size, so periodically measurements may be required again to prevent any leakage.

**Prepare your skin**
Support the skin and remove the flange carefully, do not tear off but push the skin back from the flange, use a wet gauze if necessary. Wash the skin surrounding your stoma with warm and unperfumed water. Do not use anything that contains alcohol, perfume or ether to clean your skin as it can cause damage. Do not rub the skin too hard as this might lead to skin damage. Dry the skin carefully and thoroughly to ensure that your pouch adheres well. Optional: If your skin is irritated and sore you can use skin barrier cream, hydrocolloid powder or skin barrier wipes. Apply the powder, skin barrier cream or wipe gently to the surrounding skin.

**Close the outlet**
First make sure that the closure is closed safely. Roll-up the tail of the pouch until the velcro is at the same height as the folding tabs. Press the velcro on the folding tabs firmly onto the velcro of the closure.

**Fit your pouch**
Remove the release paper of the flange. Fold the flange halfway and place the opening at the underside of the stoma. Stretch the skin of your abdomen with a flat hand and pull the whole flange over the stoma into the right position. Work your way around the stoma so that the hydrocolloid will attach to the skin and make sure that the flange is secure. Place your hand over the flange and make sure that a good seal is created.

**Empty your pouch**
Lift up the tail of the pouch and open the bag. Empty the contents of your pouch into the toilet.

**Vent**
If the pouch is full of gasses (ballooning) or if there is no air in the pouch (a vacuum) the vent can be used to let air in or out of the pouch. Make sure that the pouch is not full of faeces when opening the vent. If the pouch is full than empty the pouch or change it for a fresh one. The vent can be opened by taking off the lid. Carefully close the lid after venting.

**Viewing window**
The ostomy pouch includes a viewing window which makes it possible to inspect the stoma as well as the output from the stoma. It is also very helpful when applying the pouch to make sure it is located correctly over the stoma. Lift up the non-woven of the viewing window for inspection.

**Disposal of your pouch**
Remove your pouch and dispose of it in the trash. Do not flush your used pouch down the toilet because this will block your sewage system.

**Frequency of pouch change**
Change/empty your pouch as advised by your healthcare professional.

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**Oxmed International GmbH**